




LAB-CARE INSTRUMENTS & SERVICES

Mfg. of : Laboratory, Scientific, Surgical,
Industrial Heating & Cooling Equipments & Instruments

LEAK TEST APPARATUS



LC-33-SL

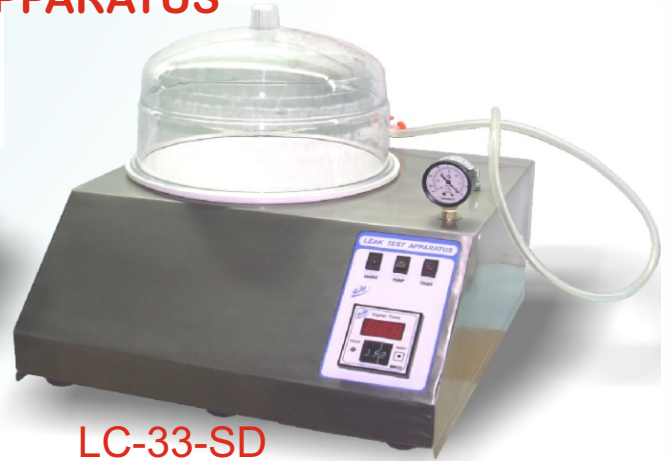

LC-33-SD

LC-33-S :Semi Automatic

Ideal instrument to check blister Packs, Strips for Leak Test. The old method of separate desiccator and huge vacuum pump is replaced with the new Leak Test Apparatus in GMP Stainless Steel Body with phenomenal compactness. The instrument is equipped with Diaphragm type vacuum pump gives higher level of vacuum in shortest possible time. Die pressed polycarbonate vacuum desiccator sustains vacuum for long time. A vacuum gauge is provided to indicate vacuum level which is connected to isolation valve to disconnect the source. Thus avoiding the need of continuous running of vacuum pump. 99 minute Digital timer helps the operator to do other work. Its amazing compactness saves a lot of valuable place in process laboratories.

UNIQUE FEATURES

- *Extremely "Compact" construction.
- *Diaphragm type vacuum pump.
- *Polycarbonate Vacuum Desiccator.
- *Stainless Steel(GMP) construction.
- *Vacuum Level upto 22" of Hg.

SPECIFICATION :

Desiccator Dia : 150/200/250/300mm

Vacuum Level : 22"Hg

Desiccator : Polycarbonate

Rating : 5A 230V single phase

Timer : 999 sec. Digital

Body : Fully stainless Steel.304

*N.B. : Due to continuous R&D model / specifications may change.
Picture shown with optional accessories.*

Specification & Dimensions are subject to modifications.